

TOSCA Laser Control System - Architecture Overview

Document Version: 1.0 **Date:** 2025-10-26 **Status:** Phase 5 - Testing & Quality Assurance **Security Status:** WARNING: **Encryption NOT Implemented** (Planned Phase 6)

Note: Current version does NOT include data encryption. See 08_security_architecture.md for planned encryption implementation (Phase 6+). DO NOT use for clinical trials or production deployment until encryption is implemented.

Executive Summary

This document outlines the architecture for a laser control system. The system integrates laser control, linear actuator positioning, GPIO-based safety interlocks, camera-based alignment, and comprehensive subject/session tracking.

System Purpose

Control and monitor laser treatments with:

- Precise power and timing control
- Adjustable ring size via linear actuator
- Real-time safety monitoring via photodiode and hotspot smoothing device
- Camera-based alignment and focus verification
- Complete treatment recording and audit trail
- Longitudinal subject tracking across multiple sessions

Technology Stack

Core Technologies

- **Language:** Python 3.10+
- **GUI Framework:** PyQt6 (modern, cross-platform, feature-rich)
- **OS Platform:** Windows 10 (Mini PC)
- **Database:** SQLite (local, single-user)

Key Libraries

```
# UI & Visualization
PyQt6           # Main GUI framework
pyqtgraph       # Real-time plotting (photodiode, power graphs)

# Image Processing & Computer Vision
opencv-python  (cv2)    # Ring detection, focus measurement
```

```

numpy                # Image array operations
pillow               # Image saving/conversion


# Hardware Interfaces
pyserial             # Arroyo laser serial communication + Arduino
                    Nano COM4
# pyfirmata           # REMOVED - replaced by custom serial
                    protocol (Oct 2025)
# Xeryon library       # Linear actuator control (existing)
# VmbPy SDK            # Allied Vision camera interface (existing)

# Database & Data Management
sqlite3              # Built-in database
sqlalchemy            # ORM for cleaner database code
alembic               # Database migrations

# Logging & Utilities
logging              # Event logging
python-dateutil       # Timestamp handling
pydantic              # Configuration validation
jsonschema           # Protocol validation
``text

```

High-Level Architecture

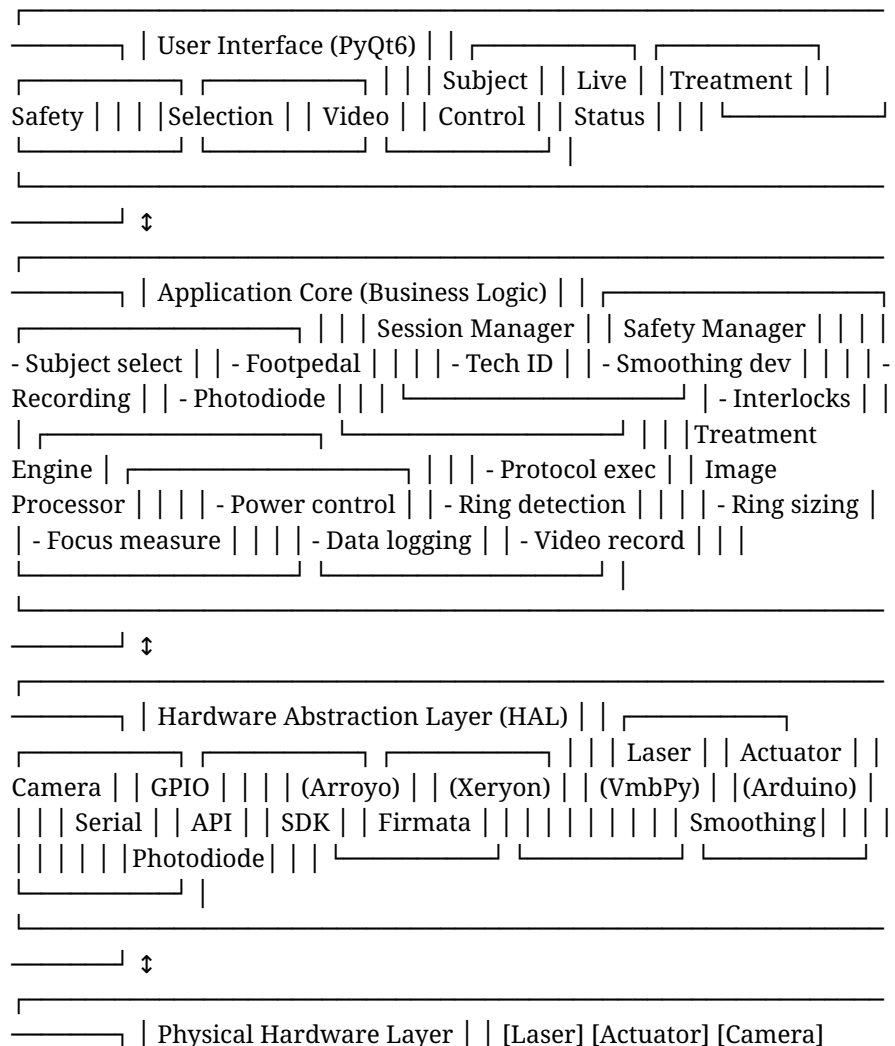

[Footpedal] [Smoothing] [Photodiode] |

Hardware Components

1. Laser Controller

- **Device:** Arroyo Instruments TEC Controller
- **Interface:** Serial communication (RS-232/USB)
- **Library:** Custom Python **class** for Arroyo serial protocol
- **Control:** Power settings, on/off, status queries

2. Linear Actuator

- **Device:** Xeryon linear stage
- **Interface:** External API library
- **Function:** Controls laser ring size by adjusting optical position
- **Control:** Position commands → Ring diameter mapping

3. Camera System

- **SDK:** VmbPy (Allied Vision Vimba Python SDK)
- **Interface:** USB/GigE
- **Functions:**
 - Live video feed display
 - Manual focus/alignment by operator
 - Ring detection (circle finding)
 - Focus quality measurement
 - Treatment recording

4. GPIO Controller - Safety Interlocks and Monitoring (Arduino Nano)

- **Device:** Arduino Nano (ATmega328P) on COM4
- **Migration Note:** Replaced FT232H GPIO expander **in** October 2025
- **Firmware:** Custom watchdog firmware **with** serial protocol
- **Communication:** USB serial (pyserial, 115200 baud)
- **Digital I/O:**
 - **Pin D2 (Output):** Smoothing device motor control
 - **Pin D3 (Input):** Smoothing device vibration sensor
 - **I2C Bus (A4/A5):** MCP4725 DAC control **for** SEMINEX aiming beam (via LDD200 driver)
- **Analog Input:**
 - **Pin A0 (ADC):** Photodiode voltage monitoring (0-5V, 10-bit)
- **Functions:**
 - Safety interlock monitoring (motor + vibration detection)
 - Real-time laser power measurement via photodiode
 - SEMINEX aiming beam control **for** alignment (12-bit DAC, 0-4095)
 - Hardware watchdog timer (1000ms timeout)
 - Cross-platform support (Windows/Linux/macOS)

Safety Architecture

Critical Safety Interlocks (All Must Pass for Laser Operation)

1. **Footpedal Deadman Switch** (GPIO-1)
 - Type: Active-high requirement
 - Behavior: Laser can only fire **while** footpedal **is** DEPRESSED
 - Fail-safe: Releasing pedal immediately disables laser
 - Poll rate: 100Hz minimum
2. **Hotspot Smoothing Device** (GPIO-1)

- Type: Signal health monitoring
 - Behavior: Device must output valid signal
 - Fail-safe: Loss of signal triggers immediate laser shutdown
 - Validation: Signal presence + value within acceptable **range**
3. ****Photodiode Feedback**** (GPIO-2 ADC)
- Type: Output power verification
 - Behavior: Measured power must match commanded power
 - Fail-safe: Deviation beyond threshold triggers shutdown
 - Monitoring: Continuous during treatment
4. ****Software E-stop****
- Type: UI button + keyboard shortcut (e.g., ESC key)
 - Behavior: Immediate treatment halt
 - Priority: Highest - bypasses **all** queues
5. ****Session Active****
- Type: Logical interlock
 - Behavior: Laser cannot fire outside active treatment session
 - Purpose: Ensures **all** actions are logged **and** attributed
6. ****Image Valid****
- Type: Camera feed health check
 - Behavior: Valid image frame received within timeout
 - Purpose: Ensures alignment/monitoring capability

Safety State Machine

[SYSTEM_OFF] → [INITIALIZING] → [READY] ↓ [FAULT] ← ← ← ← ← ← ←
 [ARMED] ← ← (all interlocks pass) ↓ ↓ [SAFE_SHUTDOWN] [TREATING] ←
 (footpedal depressed) ↓ [TREATMENT_COMPLETE]

Any interlock failure → Immediate transition to FAULT state → Safe shutdown

Session Workflow

Session Initialization

1. Application Launch ↓
2. Hardware Connection & Self-Test ↓
3. Tech ID Entry (required for all operations) ↓
4. Subject Selection Screen |— Option A: Select Existing Subject (search by subject code) | |— Load subject history |— Option B: Create New Subject |— Generate subject code, enter demographics ↓
5. Session Creation |— Log: Subject ID, Tech ID, Start Time

Pre-Treatment Setup

1. Display Live Camera Feed ↓
2. Operator Manual Actions (outside software control): |— Adjust focus (physical optics) |— Align laser ring to treatment site |— Position subject ↓
3. Software Assistance: |— Real-time focus quality indicator |— Ring detection overlay |— Alignment guides ↓
4. Operator confirms ready

Treatment Execution

1. Select Treatment Protocol |— Load saved protocol, OR |— Create/modify custom protocol ↓
2. Safety Pre-checks |— All hardware connected |— Interlocks in valid state |— Camera image valid |— Session active ↓
3. Operator initiates FIRE trigger ↓
4. System transitions to ARMED state ↓
5. Treatment Loop (while footpedal depressed): |— Execute protocol step (power, ring size) |— Monitor photodiode |— Monitor smoothing device |— Capture camera frames |— Log all parameters (timestamp, power, position, voltage) |— Check safety interlocks (every cycle) ↓
6. Treatment completion or pedal release ↓
7. Return to READY state

Session Recording

****Continuous Recording During Treatment:****

- Video: Full treatment video saved to session folder
- Event log: Every parameter change, every cycle
- Images: Periodic snapshots + key events
- Metadata: Timestamps, device states, operator actions

****Data Storage Location:****

data/ |— sessions/ | |— session_YYYYMMDD_HHMMSS/ | | |— video.avi
| | |— events.json | | |— snapshots/ | | | |— frame_001.png | | | |—
frame_002.png | | |— metadata.json

Session Closure

1. Operator ends treatment ↓
2. Save final recordings ↓
3. Add session notes ↓
4. Mark session as complete in database ↓
5. Update subject last_modified timestamp ↓
6. Return to Subject Selection (for next subject)

Treatment Protocol Engine

> ****Note:**** This section describes the older step-based protocol model.
>
> ****Current Implementation:**** See `06_protocol_builder.md` for the action-based protocol engine (current design).
>
> The action-based model provides greater flexibility with event-driven actions, conditional logic, and real-time adjustments. The step-based model shown below is kept for historical context.

Protocol Structure (Legacy Step-Based Model)

```python

```
{
 "protocol_name": "Standard Treatment A",
 "description": "5W constant for 60s at 3mm ring",
 "steps": [
 {
 "step_number": 1,
 "duration_seconds": 60,
 "power_start_watts": 5.0,
 "power_end_watts": 5.0, # Same as start = constant
 "ring_size_mm": 3.0,
 "ramp_type": "constant" # or "linear", "logarithmic"
 }
]
}
```text

```

Protocol Types

1. ****Constant Power****
 - Fixed power for duration
 - Example: 5W for 60 seconds
2. ****Linear Ramp****
 - Power increases/decreases linearly
 - Example: Ramp from 2W to 6W over 90 seconds
3. ****Multi-Step****
 - Multiple sequential steps
 - Example: 3W for 30s, then 5W for 30s, then 3W for 30s
4. ****Custom****
 - In-app protocol builder
 - Adjust any parameter on the fly

Real-Time Protocol Adjustment

- Operator can pause treatment
- Modify power/ring size during pause
- Changes logged as protocol deviation
- Resume with modified parameters

Image Processing Pipeline

Pipeline Overview

Camera Frame ↓ [Preprocessing] ─ Grayscale conversion ─ Noise reduction ─ Contrast enhancement ↓ [Ring Detection] ─ Hough Circle Transform ─ Edge detection refinement ─ Circle parameters (center, radius) ↓ [Focus Measurement] ─ Laplacian variance (sharpness) ─ Gradient magnitude ─ Focus score (0-100) ↓ [Display Overlay] ─ Detected ring outline ─ Focus indicator ─ Alignment guides ↓ [Recording] ─ Save annotated frames

Ring Detection Algorithm

```
```python
def detect_laser_ring(frame):
 """

```

```

Detect circular laser ring in camera frame

Returns:
 center (x, y), radius, confidence
 """
1. Preprocess
gray = cv2.cvtColor(frame, cv2.COLOR_BGR2GRAY)
blurred = cv2.GaussianBlur(gray, (9, 9), 2)

2. Detect circles (Hough Transform)
circles = cv2.HoughCircles(
 blurred,
 cv2.HOUGH_GRADIENT,
 dp=1,
 minDist=100,
 param1=50,
 param2=30,
 minRadius=20,
 maxRadius=200
)

3. Select best circle (brightest, most circular)
... validation logic ...

 return center, radius, confidence
```text

#### Focus Quality Measurement

```python
def calculate_focus_score(frame):
 """
 Calculate image sharpness/focus quality

 Returns:
 focus_score (0-100, higher = better focus)
 """
 gray = cv2.cvtColor(frame, cv2.COLOR_BGR2GRAY)

 # Laplacian variance method
 laplacian = cv2.Laplacian(gray, cv2.CV_64F)
 variance = laplacian.var()

 # Normalize to 0-100 scale
 focus_score = min(100, variance / 10) # Calibrate threshold

 return focus_score
```text

#### Data Architecture

#### Database: SQLite

**Location:** `data/laser_control.db`

**Key Tables:**
1. `subjects` - Subject records (anonymized)
2. `sessions` - Treatment sessions
3. `treatment_events` - Detailed event log (high frequency)
4. `protocols` - Saved treatment protocols

```

5. `calibrations` - Device calibration data
6. `safety_log` - Safety events and faults
7. `tech_users` - Technician/operator accounts

See: ``02_database_schema.md`` for full schema

Event Logging Strategy

****Two-tier logging:****

1. **High-frequency data** (100Hz+): JSON files in session folder
 - Photodiode readings
 - Camera frame metadata
 - Real-time interlock states
2. **Event-based data**: SQLite database
 - Protocol steps
 - Power changes
 - Ring size adjustments
 - Safety triggers
 - User actions

Project Directory Structure

```

laser-control-system/ ├── src/ ├── main.py # Application entry point | | |
├── config/ | | ├── settings.py # User-configurable settings | | ├──
safety_limits.py # Hard-coded safety parameters | | └──
hardware_config.py # Hardware connection parameters | | | ├── ui/ | | |
├── main_window.py # Main application window | | ├──
subject_selection.py # Subject selection/creation dialog | | ├──
treatment_control.py # Treatment control panel | | ├── video_display.py #
Live camera feed widget | | ├── protocol_builder.py # Protocol
creation/editing UI | | ├── safety_panel.py # Safety status indicators | |
└── widgets/ # Reusable UI components | | | ├── core/ | | | ├──
session_manager.py # Session lifecycle management | | ├──
treatment_engine.py # Protocol execution engine | | ├──
safety_manager.py # Safety interlock orchestration | | ├──
recording_manager.py # Video/data recording | | └──
calibration_manager.py # Calibration routines | | | ├── hardware/ | | | ├──
base.py # Abstract hardware device class | | ├── laser_controller.py #
Arroyo laser interface | | ├── actuator_controller.py # Xeryon actuator
interface | | ├── camera_controller.py # VmbPy camera interface | | ├──
gpio_interlocks.py # GPIO-1: Footpedal + Smoothing | | ├──
gpio_photodiode.py # GPIO-2: Photodiode ADC | | └──
hardware_manager.py # Unified hardware coordination | | | ├──
image_processing/ | | ├── ring_detector.py # Laser ring circle detection | | |
├── focus_analyzer.py # Focus quality measurement | | ├──
video_recorder.py # Video file writing | | └── frame_processor.py # Image
preprocessing pipeline | | | ├── database/ | | | ├── models.py #
SQLAlchemy ORM models | | ├── db_manager.py # Database operations | |
| ├── session_logger.py # High-frequency session logging | | └──
migrations/ # Alembic migration scripts | | | ├── utils/ | | | ├── logger.py #
Application logging setup | | ├── validators.py # Input validation functions | |
├── exceptions.py # Custom exception classes | └── constants.py # System-

```


wide constants | ├── data/ | ├── laser_control.db # SQLite database | ├──
sessions/ # Per-session data folders | | ├── session_/ | | ├── video.avi | |
| ├── events.json | | ├── photodiode_log.csv | | ├── snapshots/ | | ├──
metadata.json | ├── logs/ # Application logs | ├── app_YYYYMMDD.log |
| ├── errors_YYYYMMDD.log | ├── docs/ | ├── architecture/ | | ├──
01_system_overview.md # This file | | ├── 02_database_schema.md | | ├──
03_safety_system.md | | ├── 04_treatment_protocols.md | | ├──
05_image_processing.md | ├── user_manual.md | ├── installation.md |
| ├── tests/ | ├── test_hardware/ | ├── test_core/ | ├── test_safety/ | ├──
test_integration/ | ├── requirements.txt | ├── setup.py | ├── README.md ``

Development Phases

Phase 1: Foundation (Hardware + Safety)

- ☐ Hardware abstraction layer for all devices
- ☐ Safety interlock system
- ☐ Basic GUI shell (PyQt6)
- ☐ Database schema and basic CRUD operations

Phase 2: Core Treatment Features

- ☐ Subject selection and session management
- ☐ Treatment protocol engine
- ☐ Manual treatment control (constant power)
- ☐ Basic event logging

Phase 3: Advanced Features

- ☐ Ring detection and focus measurement
- ☐ Video recording
- ☐ Protocol builder UI
- ☐ Advanced ramping protocols

Phase 4: Polish & Validation

- ☐ Comprehensive testing
- ☐ User manual
- ☐ Calibration procedures
- ☐ Performance optimization

Key Design Principles

1. **Safety First:** Multiple redundant interlocks, fail-safe design
2. **Audit Trail:** Every action logged, immutable records
3. **User Workflow:** Match clinical workflow, minimize clicks
4. **Hardware Abstraction:** Easy to swap/upgrade devices
5. **Modularity:** Loosely coupled components

6. **Testability:** Unit tests for critical paths
7. **Documentation:** Code + user docs maintained together

Next Documentation Files

1. 02_database_schema.md - Complete SQL schema with indexes and constraints
2. 03_safety_system.md - Detailed safety architecture and fault handling
3. 04_treatment_protocols.md - Protocol format, execution engine, validation
4. 05_image_processing.md - Computer vision algorithms and calibration

Document Owner: System Architect **Last Updated:** 2025-10-26 **Review Frequency:** Weekly during development