

TOSCA Laser Control System - Architecture Overview

Document Version: 1.0 **Date:** 2025-10-26 **Status:** Phase 5
- Testing & Quality Assurance **Security Status:** WARNING:
Encryption NOT Implemented (Planned Phase 6)

Note: Current version does NOT include data encryption. See 08_security_architecture.md for planned encryption implementation (Phase 6+). DO NOT use for clinical trials or production deployment until encryption is implemented.

Executive Summary

This document outlines the architecture for a laser control system. The system integrates laser control, linear actuator positioning, GPIO-based safety interlocks, camera-based alignment, and comprehensive subject/session tracking.

System Purpose

Control and monitor laser treatments with:

- Precise power and timing control
- Adjustable ring size via linear actuator
- Real-time safety monitoring via photodiode and hotspot smoothing device
- Camera-based alignment and focus verification
- Complete treatment recording and audit trail
- Longitudinal subject tracking across multiple sessions

Technology Stack

Core Technologies

- **Language:** Python 3.10+
- **GUI Framework:** PyQt6 (modern, cross-platform, feature-rich)
- **OS Platform:** Windows 10 (Mini PC)
- **Database:** SQLite (local, single-user)

Key Libraries

UI & Visualization

PyQt6

Main GUI framework

```

pyqtgraph                                # Real-time plotting (photodiode, power graphs)

# Image Processing & Computer Vision
opencv-python (cv2)                      # Ring detection, focus measurement
numpy                                    # Image array operations
pillow                                   # Image saving/conversion


# Hardware Interfaces
pyserial                                # Arroyo laser serial communication + Arduino Nano COM4
# pyfirmata                             # REMOVED - replaced by custom serial protocol (Oct 2025)
# Xeryon library                         # Linear actuator control (existing)
# VmbPy SDK                             # Allied Vision camera interface (existing)


# Database & Data Management
sqlite3                                  # Built-in database
sqlalchemy                              # ORM for cleaner database code
alembic                                 # Database migrations

# Logging & Utilities
logging                                  # Event logging
python-dateutil                          # Timestamp handling
pydantic                                 # Configuration validation
jsonschema                              # Protocol validation
```text

```

## ## High-Level Architecture

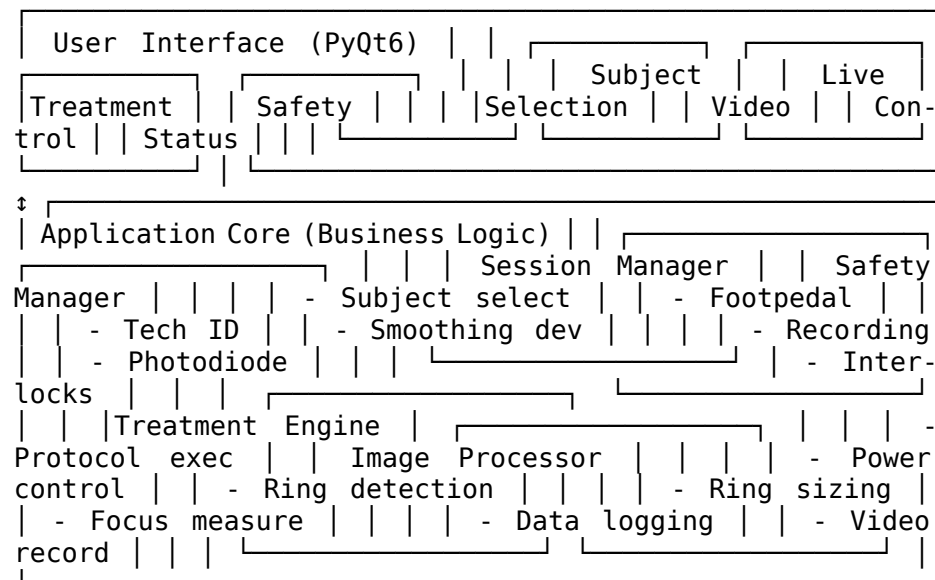

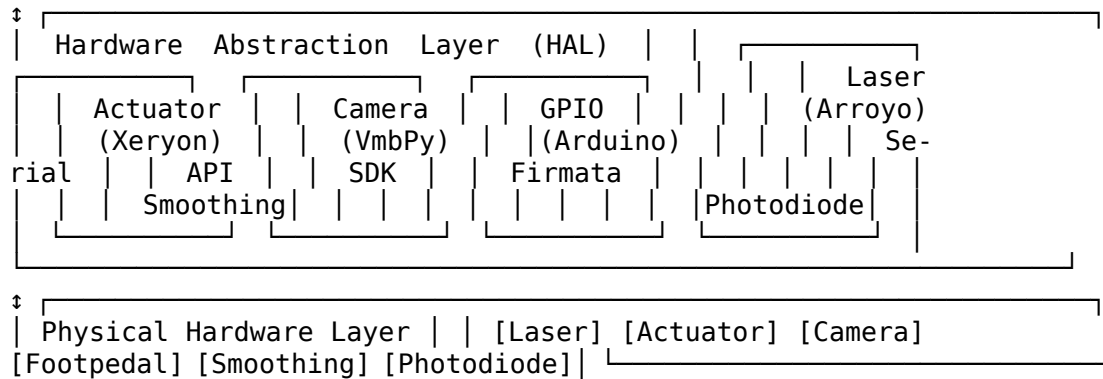

## ## Hardware Components

### ### 1. Laser Controller

- **Device:** Arroyo Instruments TEC Controller
- **Interface:** Serial communication (RS-232/USB)
- **Library:** Custom Python **class** for Arroyo serial protocol
- **Control:** Power settings, on/off, status queries

### ### 2. Linear Actuator

- **Device:** Xeryon linear stage
- **Interface:** External API library
- **Function:** Controls laser ring size by adjusting optical position
- **Control:** Position commands → Ring diameter mapping

### ### 3. Camera System

- **SDK:** VmbPy (Allied Vision Vimba Python SDK)
- **Interface:** USB/GigE
- **Functions:**
  - Live video feed display
  - Manual focus/alignment by operator
  - Ring detection (circle finding)
  - Focus quality measurement
  - Treatment recording

### ### 4. GPIO Controller - Safety Interlocks and Monitoring (Arduino Nano)

- **Device:** Arduino Nano (ATmega328P) on COM4
- **Migration Note:** Replaced FT232H GPIO expander **in** October 2025
- **Firmware:** Custom watchdog firmware **with** serial protocol
- **Communication:** USB serial (pyserial, 115200 baud)
- **Digital I/O:**
  - **Pin D2 (Output):** Smoothing device motor control
  - **Pin D3 (Input):** Smoothing device vibration sensor

- **\*\*I2C Bus (A4/A5):\*\*** MCP4725 DAC control **for** SEMINEX aiming beam (via LDD200 driver)
- **\*\*Analog Input:\*\***
  - **\*\*Pin A0 (ADC):\*\*** Photodiode voltage monitoring (**0-5V**, **10-bit**)
- **\*\*Functions:\*\***
  - Safety interlock monitoring (motor + vibration detection)
  - Real-time laser power measurement via photodiode
  - SEMINEX aiming beam control **for** alignment (**12-bit** DAC, **0-4095**)
  - Hardware watchdog timer (**1000ms** timeout)
  - Cross-platform support (Windows/Linux/macOS)

## *## Safety Architecture*

### *### Critical Safety Interlocks (All Must Pass for Laser Operation)*

1. **\*\*Footpedal Deadman Switch\*\*** (GPIO-1)
  - Type: Active-high requirement
  - Behavior: Laser can only fire **while** footpedal **is** DEPRESSED
  - Fail-safe: Releasing pedal immediately disables laser
  - Poll rate: **100Hz** minimum
2. **\*\*Hotspot Smoothing Device\*\*** (GPIO-1)
  - Type: Signal health monitoring
  - Behavior: Device must output valid signal
  - Fail-safe: Loss of signal triggers immediate laser shutdown
  - Validation: Signal presence + value within acceptable **range**
3. **\*\*Photodiode Feedback\*\*** (GPIO-2 ADC)
  - Type: Output power verification
  - Behavior: Measured power must match commanded power
  - Fail-safe: Deviation beyond threshold triggers shutdown
  - Monitoring: Continuous during treatment
4. **\*\*Software E-stop\*\***
  - Type: UI button + keyboard shortcut (e.g., ESC key)
  - Behavior: Immediate treatment halt
  - Priority: Highest - bypasses **all** queues
5. **\*\*Session Active\*\***
  - Type: Logical interlock
  - Behavior: Laser cannot fire outside active treatment session
  - Purpose: Ensures **all** actions are logged **and** attributed
6. **\*\*Image Valid\*\***
  - Type: Camera feed health check
  - Behavior: Valid image frame received within timeout
  - Purpose: Ensures alignment/monitoring capability

### ### Safety State Machine

[SYSTEM\_OFF] → [INITIALIZING] → [READY] ↓ [FAULT] ← ← ← ← ←  
← ← ← [ARMED] ← ← (all interlocks pass) ↓ ↓ [SAFE\_SHUTDOWN]  
[TREATING] ← (footpedal depressed) ↓ [TREATMENT\_COMPLETE]

Any interlock failure → Immediate transition to FAULT state → Safe shutdown

## ## Session Workflow

### ### Session Initialization

1. Application Launch ↓
2. Hardware Connection & Self-Test ↓
3. Tech ID Entry (required for all operations) ↓
4. Subject Selection Screen | Option A: Select Existing Subject (search by subject code) | | Load subject history | Option B: Create New Subject | Generate subject code, enter demographics ↓
5. Session Creation | Log: Subject ID, Tech ID, Start Time

### ### Pre-Treatment Setup

1. Display Live Camera Feed ↓
2. Operator Manual Actions (outside software control): | Adjust focus (physical optics) | Align laser ring to treatment site | Position subject ↓
3. Software Assistance: | Real-time focus quality indicator | Ring detection overlay | Alignment guides ↓
4. Operator confirms ready

### ### Treatment Execution

1. Select Treatment Protocol | Load saved protocol, OR | Create/modify custom protocol ↓
2. Safety Pre-checks | All hardware connected | Interlocks in valid state | Camera image valid | Session active ↓
3. Operator initiates FIRE trigger ↓
4. System transitions to ARMED state ↓
5. Treatment Loop (while footpedal depressed): | Execute protocol step (power, ring size) | Monitor photodiode | Monitor smoothing device | Capture camera frames |

```

 Log all parameters (timestamp, power, position, voltage)
 └─ Check safety interlocks (every cycle) ↓
6. Treatment completion or pedal release ↓
7. Return to READY state

```

### ### Session Recording

#### \*\*Continuous Recording During Treatment:\*\*

- Video: Full treatment video saved to session folder
- Event log: Every parameter change, every cycle
- Images: Periodic snapshots + key events
- Metadata: Timestamps, device states, operator actions

#### \*\*Data Storage Location:\*\*

```

data/ └─ sessions/ └─ session_YYYYMMDD_HHMMSS / └─
video.avi └─ events.json └─ snapshots/ └─
frame_001.png └─ frame_002.png └─ metadata.json

```

### ### Session Closure

1. Operator ends treatment ↓
2. Save final recordings ↓
3. Add session notes ↓
4. Mark session as complete in database ↓
5. Update subject last\_modified timestamp ↓
6. Return to Subject Selection (for next subject)

### ## Treatment Protocol Engine

> **\*\*Note:\*\*** This section describes the older step-based protocol model.

>

> **\*\*Current Implementation:\*\*** See `06\_protocol\_builder.md` for the action-based protocol engine

>

> The action-based model provides greater flexibility with event-driven actions, conditional

### ### Protocol Structure (Legacy Step-Based Model)

```
```python
```

```

{
    "protocol_name": "Standard Treatment A",
    "description": "5W constant for 60s at 3mm ring",
    "steps": [
        {
            "step_number": 1,

```

```

        "duration_seconds": 60,
        "power_start_watts": 5.0,
        "power_end_watts": 5.0, # Same as start = constant
        "ring_size_mm": 3.0,
        "ramp_type": "constant" # or "linear", "logarithmic"
    }
]
}
```text

```

### ### Protocol Types

1. **\*\*Constant Power\*\***
  - Fixed power for duration
  - Example: 5W for 60 seconds
2. **\*\*Linear Ramp\*\***
  - Power increases/decreases linearly
  - Example: Ramp from 2W to 6W over 90 seconds
3. **\*\*Multi-Step\*\***
  - Multiple sequential steps
  - Example: 3W for 30s, then 5W for 30s, then 3W for 30s
4. **\*\*Custom\*\***
  - In-app protocol builder
  - Adjust any parameter on the fly

### ### Real-Time Protocol Adjustment

- Operator can pause treatment
- Modify power/ring size during pause
- Changes logged as protocol deviation
- Resume with modified parameters

## ## Image Processing Pipeline

### ### Pipeline Overview

```

Camera Frame ↓ [Preprocessing] ─ Grayscale conversion ─
Noise reduction ─ Contrast enhancement ↓ [Ring Detection]
─ Hough Circle Transform ─ Edge detection refinement ─
Circle parameters (center, radius) ↓ [Focus Measurement]
─ Laplacian variance (sharpness) ─ Gradient magnitude ─
Focus score (0-100) ↓ [Display Overlay] ─ Detected ring
outline ─ Focus indicator ─ Alignment guides ↓ [Record-

```

ing] └ Save annotated frames

### Ring Detection Algorithm

```
```python
def detect_laser_ring(frame):
    """
    Detect circular laser ring in camera frame

    Returns:
        center (x, y), radius, confidence
    """
    # 1. Preprocess
    gray = cv2.cvtColor(frame, cv2.COLOR_BGR2GRAY)
    blurred = cv2.GaussianBlur(gray, (9, 9), 2)

    # 2. Detect circles (Hough Transform)
    circles = cv2.HoughCircles(
        blurred,
        cv2.HOUGH_GRADIENT,
        dp=1,
        minDist=100,
        param1=50,
        param2=30,
        minRadius=20,
        maxRadius=200
    )

    # 3. Select best circle (brightest, most circular)
    # ... validation logic ...

    return center, radius, confidence
```text
```

### Focus Quality Measurement

```
```python
def calculate_focus_score(frame):
    """
    Calculate image sharpness/focus quality

    Returns:
        focus_score (0-100, higher = better focus)
    """
    gray = cv2.cvtColor(frame, cv2.COLOR_BGR2GRAY)
```



```

    # Laplacian variance method
    laplacian = cv2.Laplacian(gray, cv2.CV_64F)
    variance = laplacian.var()

    # Normalize to 0-100 scale
    focus_score = min(100, variance / 10) # Calibrate threshold

    return focus_score
```text

Data Architecture

Database: SQLite

Location: `data/laser_control.db`

Key Tables:
1. `subjects` - Subject records (anonymized)
2. `sessions` - Treatment sessions
3. `treatment_events` - Detailed event log (high frequency)
4. `protocols` - Saved treatment protocols
5. `calibrations` - Device calibration data
6. `safety_log` - Safety events and faults
7. `tech_users` - Technician/operator accounts

See: `02_database_schema.md` for full schema

Event Logging Strategy

Two-tier logging:

1. **High-frequency data** (100Hz+): JSON files in session folder
 - Photodiode readings
 - Camera frame metadata
 - Real-time interlock states

2. **Event-based data**: SQLite database
 - Protocol steps
 - Power changes
 - Ring size adjustments
 - Safety triggers
 - User actions

Project Directory Structure

```



README.md ““

## Development Phases

### Phase 1: Foundation (Hardware + Safety)

- ☐ Hardware abstraction layer for all devices
- ☐ Safety interlock system
- ☐ Basic GUI shell (PyQt6)
- ☐ Database schema and basic CRUD operations

### Phase 2: Core Treatment Features

- ☐ Subject selection and session management
- ☐ Treatment protocol engine
- ☐ Manual treatment control (constant power)
- ☐ Basic event logging

### Phase 3: Advanced Features

- ☐ Ring detection and focus measurement
- ☐ Video recording
- ☐ Protocol builder UI
- ☐ Advanced ramping protocols

### Phase 4: Polish & Validation

- ☐ Comprehensive testing
- ☐ User manual
- ☐ Calibration procedures
- ☐ Performance optimization

## Key Design Principles

1. **Safety First:** Multiple redundant interlocks, fail-safe design
2. **Audit Trail:** Every action logged, immutable records
3. **User Workflow:** Match clinical workflow, minimize clicks
4. **Hardware Abstraction:** Easy to swap/upgrade devices
5. **Modularity:** Loosely coupled components
6. **Testability:** Unit tests for critical paths
7. **Documentation:** Code + user docs maintained together

## **Next Documentation Files**

1. 02\_database\_schema.md - Complete SQL schema with indexes and constraints
2. 03\_safety\_system.md - Detailed safety architecture and fault handling
3. 04\_treatment\_protocols.md - Protocol format, execution engine, validation
4. 05\_image\_processing.md - Computer vision algorithms and calibration

---

**Document Owner:** System Architect **Last Updated:** 2025-10-26  
**Review Frequency:** Weekly during development