

TOSCA Laser Control System - Architecture Overview

Document Version: 1.0 **Date:** 2025-10-26 **Status:** Phase 5 - Testing & Quality Assurance **Security Status:** WARNING: Encryption NOT Implemented (Planned Phase 6)

Note: Current version does NOT include data encryption. See 08_security_architecture.md for planned encryption implementation (Phase 6+). DO NOT use for clinical trials or production deployment until encryption is implemented.

Executive Summary

This document outlines the architecture for a laser control system. The system integrates laser control, linear actuator positioning, GPIO-based safety interlocks, camera-based alignment, and comprehensive subject/session tracking.

System Purpose

Control and monitor laser treatments with:

- Precise power and timing control
- Adjustable ring size via linear actuator
- Real-time safety monitoring via photodiode and hotspot smoothing device
- Camera-based alignment and focus verification
- Complete treatment recording and audit trail
- Longitudinal subject tracking across multiple sessions

Technology Stack

Core Technologies

- **Language:** Python 3.10+
- **GUI Framework:** PyQt6 (modern, cross-platform, feature-rich)
- **OS Platform:** Windows 10 (Mini PC)
- **Database:** SQLite (local, single-user)

Key Libraries

UI & Visualization

PyQt6	<i># Main GUI framework</i>
pyqtgraph	<i># Real-time plotting (photodiode, power graphs)</i>

Image Processing & Computer Vision

opencv-python (cv2)	<i># Ring detection, focus measurement</i>
numpy	<i># Image array operations</i>
pillow	<i># Image saving/conversion</i>

Hardware Interfaces

pyserial	<i># Arroyo laser serial communication + Arduino Nano COM4</i>
<i># pyfirmata</i>	<i># REMOVED - replaced by custom serial protocol (Oct 2025)</i>
<i># Xeryon library</i>	<i># Linear actuator control (existing)</i>
<i># VmbPy SDK</i>	<i># Allied Vision camera interface (existing)</i>


Database & Data Management

```
sqlite3          # Built-in database
sqlalchemy       # ORM for cleaner database code
alembic          # Database migrations
```

Logging & Utilities

```
logging          # Event logging
python-dateutil  # Timestamp handling
pydantic         # Configuration validation
jsonschema       # Protocol validation
```text
```

## ## High-Level Architecture

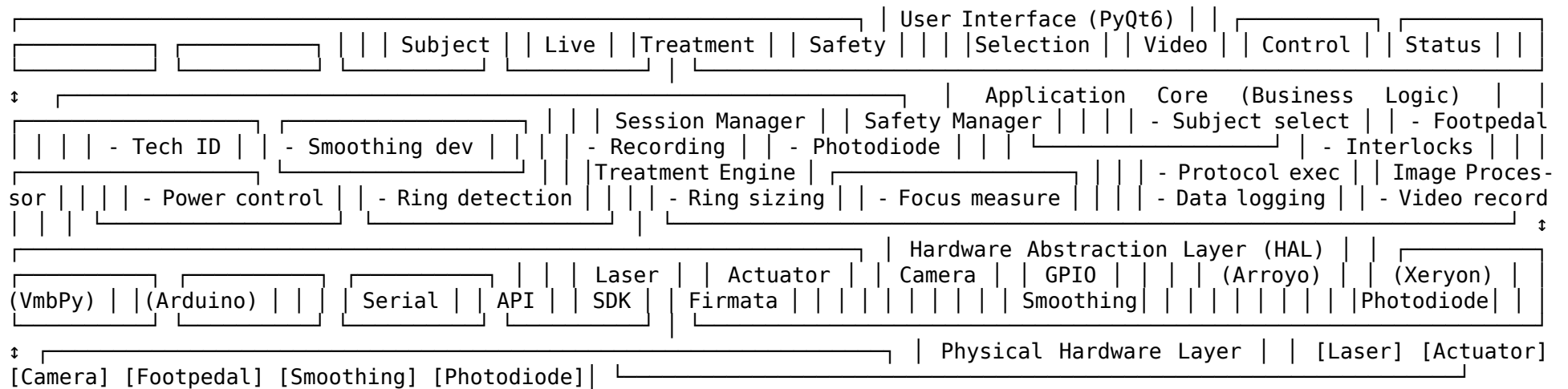

## ## Hardware Components

### ### 1. Laser Controller

- **\*\*Device:\*\*** Arroyo Instruments TEC Controller
- **\*\*Interface:\*\*** Serial communication (RS-232/USB)
- **\*\*Library:\*\*** Custom Python **class** for Arroyo serial protocol
- **\*\*Control:\*\*** Power settings, on/off, status queries

### ### 2. Linear Actuator

- **\*\*Device:\*\*** Xeryon linear stage
- **\*\*Interface:\*\*** External API library
- **\*\*Function:\*\*** Controls laser ring size by adjusting optical position
- **\*\*Control:\*\*** Position commands → Ring diameter mapping

### ### 3. Camera System

- **\*\*SDK:\*\*** VmbPy (Allied Vision Vimba Python SDK)
- **\*\*Interface:\*\*** USB/GigE
- **\*\*Functions:\*\***
  - Live video feed display
  - Manual focus/alignment by operator
  - Ring detection (circle finding)
  - Focus quality measurement
  - Treatment recording

### ### 4. GPIO Controller - Safety Interlocks and Monitoring (Arduino Nano)

- **\*\*Device:\*\*** Arduino Nano (ATmega328P) on COM4
- **\*\*Migration Note:\*\*** Replaced FT232H GPIO expander in October 2025
- **\*\*Firmware:\*\*** Custom watchdog firmware with serial protocol
- **\*\*Communication:\*\*** USB serial (pyserial, 115200 baud)
- **\*\*Digital I/O:\*\***
  - **\*\*Pin D2 (Output):\*\*** Smoothing device motor control
  - **\*\*Pin D3 (Input):\*\*** Smoothing device vibration sensor
  - **\*\*I2C Bus (A4/A5):\*\*** MCP4725 DAC control for SEMINEX aiming beam (via LDD200 driver)
- **\*\*Analog Input:\*\***
  - **\*\*Pin A0 (ADC):\*\*** Photodiode voltage monitoring (0-5V, 10-bit)
- **\*\*Functions:\*\***
  - Safety interlock monitoring (motor + vibration detection)
  - Real-time laser power measurement via photodiode
  - SEMINEX aiming beam control for alignment (12-bit DAC, 0-4095)
  - Hardware watchdog timer (1000ms timeout)
  - Cross-platform support (Windows/Linux/macOS)

## ## Safety Architecture

### ### Critical Safety Interlocks (All Must Pass for Laser Operation)

1. **\*\*Footpedal Deadman Switch\*\*** (GPIO-1)
  - Type: Active-high requirement
  - Behavior: Laser can only fire while footpedal is DEPRESSED
  - Fail-safe: Releasing pedal immediately disables laser
  - Poll rate: 100Hz minimum
2. **\*\*Hotspot Smoothing Device\*\*** (GPIO-1)
  - Type: Signal health monitoring
  - Behavior: Device must output valid signal
  - Fail-safe: Loss of signal triggers immediate laser shutdown
  - Validation: Signal presence + value within acceptable range
3. **\*\*Photodiode Feedback\*\*** (GPIO-2 ADC)

- Type: Output power verification
- Behavior: Measured power must match commanded power
- Fail-safe: Deviation beyond threshold triggers shutdown
- Monitoring: Continuous during treatment

#### 4. **\*\*Software E-stop\*\***

- Type: UI button + keyboard shortcut (e.g., ESC key)
- Behavior: Immediate treatment halt
- Priority: Highest - bypasses **all** queues

#### 5. **\*\*Session Active\*\***

- Type: Logical interlock
- Behavior: Laser cannot fire outside active treatment session
- Purpose: Ensures **all** actions are logged **and** attributed

#### 6. **\*\*Image Valid\*\***

- Type: Camera feed health check
- Behavior: Valid image frame received within timeout
- Purpose: Ensures alignment/monitoring capability

### ### *Safety State Machine*

[SYSTEM\_OFF] → [INITIALIZING] → [READY] ↓ [FAULT] ← ← ← ← ← ← ← ← [ARMED] ← ← (all interlocks pass) ↓ ↓ [SAFE\_SHUTDOWN]  
 [TREATING] ← (footpedal depressed) ↓ [TREATMENT\_COMPLETE]

Any interlock failure → Immediate transition to FAULT state → Safe shutdown

### ## Session Workflow

#### ### Session Initialization

1. Application Launch ↓
2. Hardware Connection & Self-Test ↓
3. Tech ID Entry (required for all operations) ↓
4. Subject Selection Screen | Option A: Select Existing Subject (search by subject code) | ↳ Load subject history ↳  
 Option B: Create New Subject ↳ Generate subject code, enter demographics ↓
5. Session Creation ↳ Log: Subject ID, Tech ID, Start Time

#### ### Pre-Treatment Setup

1. Display Live Camera Feed ↓
2. Operator Manual Actions (outside software control): | Adjust focus (physical optics) | Align laser ring to treatment site ↳ Position subject ↓
3. Software Assistance: | Real-time focus quality indicator | Ring detection overlay ↳ Alignment guides ↓
4. Operator confirms ready

### ### Treatment Execution

1. Select Treatment Protocol | Load saved protocol, OR | Create/modify custom protocol ↓
2. Safety Pre-checks | All hardware connected | Interlocks in valid state | Camera image valid | Session active ↓
3. Operator initiates FIRE trigger ↓
4. System transitions to ARMED state ↓
5. Treatment Loop (while footpedal depressed): | Execute protocol step (power, ring size) | Monitor photodiode | Monitor smoothing device | Capture camera frames | Log all parameters (timestamp, power, position, voltage) | Check safety interlocks (every cycle) ↓
6. Treatment completion or pedal release ↓
7. Return to READY state

### ### Session Recording

#### \*\*Continuous Recording During Treatment:\*\*

- Video: Full treatment video saved to session folder
- Event log: Every parameter change, every cycle
- Images: Periodic snapshots + key events
- Metadata: Timestamps, device states, operator actions

#### \*\*Data Storage Location:\*\*

data/ | sessions/ | | session\_YYYYMMDD\_HHMMSS/ | | | video.avi | | | events.json | | | snapshots/ | | | | frame\_001.png | | | | frame\_002.png | | | metadata.json

### ### Session Closure

1. Operator ends treatment ↓
2. Save final recordings ↓
3. Add session notes ↓
4. Mark session as complete in database ↓
5. Update subject last\_modified timestamp ↓
6. Return to Subject Selection (for next subject)

### ## Treatment Protocol Engine

> \*\*Note:\*\* This section describes the older step-based protocol model.

>

> \*\*Current Implementation:\*\* See `06\_protocol\_builder.md` for the action-based protocol engine (current design).

>

> The action-based model provides greater flexibility with event-driven actions, conditional logic, and real-time adjustments. The step-based model

### ### Protocol Structure (Legacy Step-Based Model)

```

```python
{
  "protocol_name": "Standard Treatment A",
  "description": "5W constant for 60s at 3mm ring",
  "steps": [
    {
      "step_number": 1,
      "duration_seconds": 60,
      "power_start_watts": 5.0,
      "power_end_watts": 5.0, # Same as start = constant
      "ring_size_mm": 3.0,
      "ramp_type": "constant" # or "linear", "logarithmic"
    }
  ]
}
```text

```

### ### Protocol Types

1. **\*\*Constant Power\*\***
  - Fixed power for duration
  - Example: 5W for 60 seconds
2. **\*\*Linear Ramp\*\***
  - Power increases/decreases linearly
  - Example: Ramp from 2W to 6W over 90 seconds
3. **\*\*Multi-Step\*\***
  - Multiple sequential steps
  - Example: 3W for 30s, then 5W for 30s, then 3W for 30s
4. **\*\*Custom\*\***
  - In-app protocol builder
  - Adjust any parameter on the fly

### ### Real-Time Protocol Adjustment

- Operator can pause treatment
- Modify power/ring size during pause
- Changes logged as protocol deviation
- Resume with modified parameters

### ## Image Processing Pipeline

### ### Pipeline Overview

Camera Frame ↓ [Preprocessing] ┤ Grayscale conversion ┤ Noise reduction ┤ Contrast enhancement ↓ [Ring Detection]  
 ┤ Hough Circle Transform ┤ Edge detection refinement ┤ Circle parameters (center, radius) ↓ [Focus Measurement] ┤  
 Laplacian variance (sharpness) ┤ Gradient magnitude ┤ Focus score (0-100) ↓ [Display Overlay] ┤ Detected ring outline  
 ┤ Focus indicator ┤ Alignment guides ↓ [Recording] ┤ Save annotated frames

### Ring Detection Algorithm

```
```python
def detect_laser_ring(frame):
    """
    Detect circular laser ring in camera frame

    Returns:
        center (x, y), radius, confidence
    """
    # 1. Preprocess
    gray = cv2.cvtColor(frame, cv2.COLOR_BGR2GRAY)
    blurred = cv2.GaussianBlur(gray, (9, 9), 2)

    # 2. Detect circles (Hough Transform)
    circles = cv2.HoughCircles(
        blurred,
        cv2.HOUGH_GRADIENT,
        dp=1,
        minDist=100,
        param1=50,
        param2=30,
        minRadius=20,
        maxRadius=200
    )

    # 3. Select best circle (brightest, most circular)
    # ... validation logic ...

    return center, radius, confidence
```text
```

### Focus Quality Measurement

```
```python
def calculate_focus_score(frame):
    """
    Calculate image sharpness/focus quality
```

```

Returns:
    focus_score (0-100, higher = better focus)
"""
gray = cv2.cvtColor(frame, cv2.COLOR_BGR2GRAY)

# Laplacian variance method
laplacian = cv2.Laplacian(gray, cv2.CV_64F)
variance = laplacian.var()

# Normalize to 0-100 scale
focus_score = min(100, variance / 10) # Calibrate threshold

return focus_score
```text

Data Architecture

Database: SQLite

Location: `data/laser_control.db`

Key Tables:
1. `subjects` - Subject records (anonymized)
2. `sessions` - Treatment sessions
3. `treatment_events` - Detailed event log (high frequency)
4. `protocols` - Saved treatment protocols
5. `calibrations` - Device calibration data
6. `safety_log` - Safety events and faults
7. `tech_users` - Technician/operator accounts

See: `02_database_schema.md` for full schema

Event Logging Strategy

Two-tier logging:

1. **High-frequency data** (100Hz+): JSON files in session folder
 - Photodiode readings
 - Camera frame metadata
 - Real-time interlock states

2. **Event-based data**: SQLite database
 - Protocol steps
 - Power changes
 - Ring size adjustments

```



- Safety triggers
- User actions

## ## Project Directory Structure

```

laser-control-system/ ├── src/ | ├── main.py # Application entry point | | | ├── config/ | | | ├── settings.py # User-
configurable settings | | | ├── safety_limits.py # Hard-coded safety parameters | | | └── hardware_config.py # Hardware
connection parameters | | | ├── ui/ | | | ├── main_window.py # Main application window | | | ├── subject_selection.py #
Subject selection/creation dialog | | | ├── treatment_control.py # Treatment control panel | | | ├── video_display.py #
Live camera feed widget | | | ├── protocol_builder.py # Protocol creation/editing UI | | | ├── safety_panel.py # Safety
status indicators | | | └── widgets/ # Reusable UI components | | | ├── core/ | | | ├── session_manager.py # Session
lifecycle management | | | ├── treatment_engine.py # Protocol execution engine | | | ├── safety_manager.py # Safety interlock
orchestration | | | ├── recording_manager.py # Video/data recording | | | └── calibration_manager.py # Calibration routines
| | | ├── hardware/ | | | ├── base.py # Abstract hardware device class | | | ├── laser_controller.py # Arroyo laser interface
| | | ├── actuator_controller.py # Xeryon actuator interface | | | ├── camera_controller.py # VmbPy camera interface | |
| ├── gpio_interlocks.py # GPIO-1: Footpedal + Smoothing | | | ├── gpio_photodiode.py # GPIO-2: Photodiode ADC | | | └──
hardware_manager.py # Unified hardware coordination | | | ├── image_processing/ | | | ├── ring_detector.py # Laser ring
circle detection | | | ├── focus_analyzer.py # Focus quality measurement | | | ├── video_recorder.py # Video file writing |
| | | ├── frame_processor.py # Image preprocessing pipeline | | | ├── database/ | | | ├── models.py # SQLAlchemy ORM models |
| | | ├── db_manager.py # Database operations | | | ├── session_logger.py # High-frequency session logging | | | └── migrations/
Alembic migration scripts | | | └── utils/ | | | ├── logger.py # Application logging setup | | | ├── validators.py # Input
validation functions | | | ├── exceptions.py # Custom exception classes | | | ├── constants.py # System-wide constants | | | ├── data/
| | | ├── laser_control.db # SQLite database | | | ├── sessions/ # Per-session data folders | | | ├── session_/ | | | ├── video.avi
| | | ├── events.json | | | ├── photodiode_log.csv | | | ├── snapshots/ | | | ├── metadata.json | | | ├── logs/ # Application logs |
| | | ├── app_YYYYMMDD.log | | | ├── errors_YYYYMMDD.log | | | ├── docs/ | | | ├── architecture/ | | | ├── 01_system_overview.md # This file |
| | | ├── 02_database_schema.md | | | ├── 03_safety_system.md | | | ├── 04_treatment_protocols.md | | | ├── 05_image_processing.md
| | | ├── user_manual.md | | | ├── installation.md | | | ├── tests/ | | | ├── test_hardware/ | | | ├── test_core/ | | | ├── test_safety/ | |
| └── test_integration/ | | | ├── requirements.txt | | | ├── setup.py | | | └── README.md ""

```

## Development Phases

### Phase 1: Foundation (Hardware + Safety)

- ☐ Hardware abstraction layer for all devices
- ☐ Safety interlock system
- ☐ Basic GUI shell (PyQt6)
- ☐ Database schema and basic CRUD operations

### Phase 2: Core Treatment Features

- ☐ Subject selection and session management
- ☐ Treatment protocol engine
- ☐ Manual treatment control (constant power)
- ☐ Basic event logging

### Phase 3: Advanced Features

- ☐ Ring detection and focus measurement
- ☐ Video recording
- ☐ Protocol builder UI
- ☐ Advanced ramping protocols

### Phase 4: Polish & Validation

- ☐ Comprehensive testing
- ☐ User manual
- ☐ Calibration procedures
- ☐ Performance optimization

### Key Design Principles

1. **Safety First:** Multiple redundant interlocks, fail-safe design
2. **Audit Trail:** Every action logged, immutable records
3. **User Workflow:** Match clinical workflow, minimize clicks
4. **Hardware Abstraction:** Easy to swap/upgrade devices
5. **Modularity:** Loosely coupled components
6. **Testability:** Unit tests for critical paths
7. **Documentation:** Code + user docs maintained together

### Next Documentation Files

1. 02\_database\_schema.md - Complete SQL schema with indexes and constraints
2. 03\_safety\_system.md - Detailed safety architecture and fault handling
3. 04\_treatment\_protocols.md - Protocol format, execution engine, validation
4. 05\_image\_processing.md - Computer vision algorithms and calibration

---

**Document Owner:** System Architect **Last Updated:** 2025-10-26 **Review Frequency:** Weekly during development