

# Recording Manager Architecture

- [Architecture Diagrams](#)
- [Table of Contents](#)
- [Overview](#)
  - [Purpose](#)
- [Recording Architecture](#)
  - [System Overview](#)
  - [Key Components](#)
- [Video Format & Codec](#)
  - [Current Implementation \(Phase 5\)](#)
  - [File Naming Convention](#)
- [Session Integration](#)
  - [Recording Workflow](#)
  - [CameraController API](#)
- [Future Encryption](#)
  - [Encryption Strategy \(Phase 6\)](#)
  - [Playback with Decryption](#)
- [Testing](#)
  - [Recording Tests](#)
- [Performance Considerations](#)
  - [Disk Space](#)
  - [CPU Usage](#)
  - [Memory Usage](#)
- [References](#)
  - [OpenCV](#)
  - [Video Encryption](#)

## Architecture Diagrams

**Document Version:** 1.0 **Last Updated:** 2025-10-26 **Status:** Implemented - Phase 5 **WARNING: Videos NOT Encrypted** (Encryption planned Phase 6) **Priority:** HIGH - Required for treatment documentation

---

**WARNING: WARNING - Videos Stored Unencrypted:**

**Current Implementation:** Video recordings are stored as **plaintext MP4 files** without encryption.

- [FAILED] Videos are **NOT encrypted** (HIPAA vulnerability)
- [FAILED] PHI/PII visible in video (patient faces, treatment areas)
- [FAILED] DO NOT use for clinical trials (encryption required)

**Future Implementation (Phase 6):** AES-256-GCM encryption (see section below)

---

## Table of Contents

1. [Overview](#)
  2. [Recording Architecture](#)
  3. [Video Format & Codec](#)
  4. [Session Integration](#)
  5. [Future Encryption](#)
-

# Overview

## Purpose


TOSCA records complete video of every treatment session for documentation, quality assurance, and regulatory compliance. This document describes the recording architecture and integration with the session management system.

### Recording Objectives:

1. **Documentation:** Visual record of treatment for patient file
  2. **Quality Assurance:** Post-treatment review by clinician
  3. **Training:** Educational resource for new operators
  4. **Compliance:** FDA requires complete treatment records
  5. **Legal Protection:** Evidence in case of adverse events
- 

# Recording Architecture

## System Overview

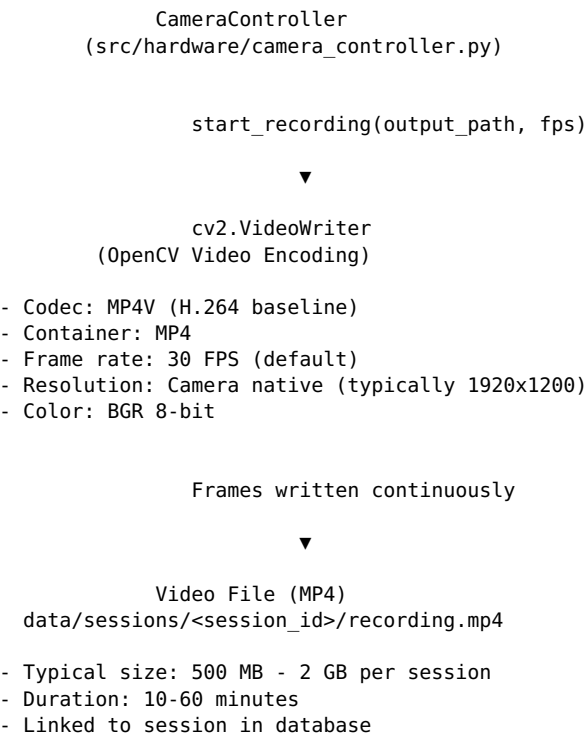

## Key Components

**CameraController** (src/hardware/camera\_controller.py) - Manages VideoWriter lifecycle - Handles frame capture and encoding - Emits signals for recording status

**Session** (src/core/session.py) - Stores video file path (video\_file\_path field) - Links recording to patient/operator - Tracks session metadata

**VideoWriter** (cv2.VideoWriter from OpenCV) - Encodes frames to MP4 - Handles codec and compression - Thread-safe writes

---

# Video Format & Codec

## Current Implementation (Phase 5)

**Container:** MP4 (.mp4)

**Video Codec:** MP4V (MPEG-4 Part 2)

```
fourcc = cv2.VideoWriter_fourcc(*"mp4v")
writer = cv2.VideoWriter(output_path, fourcc, fps, frame_size)
```

**Frame Rate:** 30 FPS (configurable)

**Resolution:** Native camera resolution (typically 1920x1200)

**Color Format:** BGR 8-bit (OpenCV default)

**Compression:** Lossy (H.264 baseline profile)

## File Naming Convention

data/sessions/<session\_id>/recording.mp4

**Example:**

data/sessions/2025-10-26\_143015\_subject-042/recording.mp4

**Metadata in Database:**

```
-- sessions table
UPDATE sessions
SET video_path = 'data/sessions/2025-10-26_143015_subject-042/recording.mp4',
    video_size_bytes = 1234567890,
    video_duration_s = 1245.5
WHERE id = 42;
```

---

# Session Integration

## Recording Workflow

### 1. Session Start:

```
# Create session
session = Session(
    session_id="2025-10-26_143015_subject-042",
    subject_id="042",
    operator_id="tech-005",
    start_time=datetime.now()
)

# Start recording
output_path = f"data/sessions/{session.session_id}/recording.mp4"
camera.start_recording(output_path, fps=30)

session.video_file_path = output_path
session.start()
```

### 2. During Treatment:

```
# Frames automatically captured and encoded
# VideoWriter writes to file continuously

# Recording status available via signal
camera.recording_status_changed.connect(update_ui)
```

### 3. Session End:

```

# Stop recording
camera.stop_recording()

# Finalize session
session.complete()

# Save to database
db.save_session(session)

```

## CameraController API

### start\_recording()

```

def start_recording(
    self,
    output_path: Path,
    fps: float = 30.0,
    codec: str = "mp4v"
) -> bool:
    """
    Start video recording to file.

    Args:
        output_path: Path to output MP4 file
        fps: Frame rate (default 30 FPS)
        codec: FourCC codec (default "mp4v" for H.264)

    Returns:
        True if recording started successfully, False otherwise

    Emits:
        recording_started signal
    """
    # Create output directory if needed
    output_path.parent.mkdir(parents=True, exist_ok=True)

    # Initialize VideoWriter
    fourcc = cv2.VideoWriter_fourcc(*codec)
    self.writer = cv2.VideoWriter(
        str(output_path),
        fourcc,
        fps,
        self.frame_size # (width, height)
    )

    if not self.writer.isOpened():
        logger.error(f"Failed to open VideoWriter: {output_path}")
        return False

    self.recording = True
    self.recording_started.emit(str(output_path))
    logger.info(f"Recording started: {output_path} @ {fps} FPS")
    return True

```

### stop\_recording()

```

def stop_recording(self) -> None:
    """
    Stop video recording and finalize file.

    Emits:
        recording_stopped signal
    """
    if not self.recording or not self.writer:
        return

    # Release VideoWriter (finalizes file)
    self.writer.release()
    self.writer = None

    self.recording = False
    self.recording_stopped.emit()
    logger.info("Recording stopped")

```

## Frame Capture (Automatic):

```
# In CameraStreamThread.run():
def frame_callback(cam, stream, frame):
    """Called for every camera frame."""
    frame_data = frame.as_numpy_ndarray()

    # If recording, write frame to video
    if self.writer and self.recording:
        self.writer.write(frame_data) # ← Thread-safe

    # Also emit for live display
    self.frame_ready.emit(frame_data)
```

---

## Future Encryption

### Encryption Strategy (Phase 6)

**See:** docs/architecture/08\_security\_architecture.md for complete details

**Summary:** - **Algorithm:** AES-256-GCM (authenticated encryption) - **Key Derivation:** PBKDF2-HMAC-SHA256 from master password - **File Format:** [12-byte nonce][encrypted video][16-byte auth tag]

**Benefits:** - [DONE] Protects PHI (patient video is sensitive data) - [DONE] Tamper-evident (authentication tag) - [DONE] HIPAA compliant (encryption at rest)

### Implementation:

```
# Phase 6 implementation
def encrypt_video_file(input_path: Path, output_path: Path, key: bytes) -> None:
    """Encrypt video file with AES-256-GCM."""
    from cryptography.hazmat.primitives.ciphers.aead import AESGCM
    import os

    aesgcm = AESGCM(key) # 32-byte key for AES-256
    nonce = os.urandom(12) # 96-bit nonce

    with open(input_path, 'rb') as f:
        plaintext = f.read()

    ciphertext = aesgcm.encrypt(nonce, plaintext, None)

    with open(output_path, 'wb') as f:
        f.write(nonce)
        f.write(ciphertext) # Includes 16-byte auth tag

# Usage during session end
camera.stop_recording()
encrypt_video_file(
    input_path=session.video_file_path,
    output_path=session.video_file_path.with_suffix('.enc'),
    key=derived_video_key
)
os.remove(session.video_file_path) # Delete plaintext
```

## Playback with Decryption

### Planned UI (Phase 6):

```
def playback_encrypted_video(encrypted_path: Path, key: bytes) -> None:
    """Decrypt and playback video."""
    from cryptography.hazmat.primitives.ciphers.aead import AESGCM

    # Read encrypted file
    with open(encrypted_path, 'rb') as f:
        nonce = f.read(12)
        ciphertext = f.read()
```

```
# Decrypt
aesgcm = AESGCM(key)
plaintext = aesgcm.decrypt(nonce, ciphertext, None)

# Write to temporary file for playback
temp_path = Path("temp/playback.mp4")
with open(temp_path, 'wb') as f:
    f.write(plaintext)

# Play using cv2.VideoCapture or external player
play_video(temp_path)

# Delete temporary file
temp_path.unlink()
```

---

## Testing

### Recording Tests

#### Test 1: Verify Recording Start/Stop

```
def test_recording_lifecycle():
    """Test complete recording workflow."""
    camera = CameraController()
    output_path = Path("test_output.mp4")

    # Start recording
    result = camera.start_recording(output_path, fps=30)
    assert result is True
    assert camera.recording is True

    # Simulate frames
    time.sleep(2.0)

    # Stop recording
    camera.stop_recording()
    assert camera.recording is False
    assert output_path.exists()

    # Verify file is valid MP4
    cap = cv2.VideoCapture(str(output_path))
    assert cap.isOpened()
    cap.release()
```

#### Test 2: Verify Frame Count

```
def test_frame_count_accuracy():
    """Verify all frames are written."""
    camera = CameraController()
    output_path = Path("test_frames.mp4")

    camera.start_recording(output_path, fps=30)
    time.sleep(3.0) # Record for 3 seconds
    camera.stop_recording()

    # Count frames in video
    cap = cv2.VideoCapture(str(output_path))
    frame_count = int(cap.get(cv2.CAP_PROP_FRAME_COUNT))
    cap.release()

    # Expect ~90 frames (30 FPS * 3 seconds)
    assert 85 <= frame_count <= 95 # Allow ±5 frames tolerance
```

---

## Performance Considerations

### Disk Space

## Estimates:

### Session Duration Resolution Frame Rate File Size

10 minutes	1920x1200	30 FPS	~500 MB
30 minutes	1920x1200	30 FPS	~1.5 GB
60 minutes	1920x1200	30 FPS	~3.0 GB

**Storage Requirements:** - Average session: 30 minutes = 1.5 GB - Daily capacity (8 sessions): ~12 GB - Weekly capacity (40 sessions): ~60 GB - Monthly capacity (160 sessions): ~240 GB

**Recommendation:** Minimum 1 TB storage for 6 months of sessions

## CPU Usage

**VideoWriter CPU Impact:** - Encoding: ~5-10% CPU (H.264 hardware acceleration if available) - Writing: ~1-2% CPU (disk I/O)

**Total:** <15% CPU overhead during recording (acceptable)

## Memory Usage

**Frame Buffers:** - Frame size: 1920x1200 x 3 bytes (BGR) = 6.9 MB - Typical buffer: 2-3 frames = ~20 MB - VideoWriter buffer: ~10 MB

**Total:** ~30 MB memory overhead (negligible)

---

## References

### OpenCV

- **VideoWriter:** [https://docs.opencv.org/4.x/dd/d43/tutorial\\_py\\_video\\_display.html](https://docs.opencv.org/4.x/dd/d43/tutorial_py_video_display.html)
- **Video Codecs:** [https://docs.opencv.org/4.x/dd/d9e/classcv\\_1\\_1VideoWriter.html](https://docs.opencv.org/4.x/dd/d9e/classcv_1_1VideoWriter.html)

### Video Encryption

- **AES-GCM:** NIST SP 800-38D
  - **See:** docs/architecture/08\_security\_architecture.md
- 

**Document Owner:** Software Architect **Last Updated:** 2025-10-26 **Next Review:** Before Phase 6 (encryption implementation) **Status:** Recording implemented - Encryption planned